

Protein information sheet



| | |
|-------------------|------------------------|
| Gene name | VCP |
| Uniprot ID | P55072 |
| Region | M1 - G806 |

| | |
|-------------------------|--|
| Description | Valosin containing protein, Transitional endoplasmic reticulum ATPase |
| Synonyms | CDC48, TERA, p97 |
| Construct ID | VCPA-c000 |
| Protein sequence | <p>SMASGADSKGDDLSTAILKQKNRPNRLIVDEAINEDNSVVSLSQPKMDELQLF RGDTVLLKGGKRREAVCIVLSDDTCSDEKIRMNRVVRNNLRVRLGDVISIQPC PDVKYGKRIHVLPIDDTVEGITGNLFEVYLKPYFLEAYRPIRKGDIFLVRGGM RAVEFKVVETDPSYCIVAPDTVIHCEGEPIKREDEEESLNEVGYYDDIGGCRK QLAQIKEMVELPLRHPALFKAIGVKPPRGILLYGPPGTGKTLIARAVANETGA FFFLINGPEIMSKLAGESESNLKAFEEAEKNAPAIIFIDELDAIAPKREKTH GEVERRIVSQLLTMDGLKQRAHVIVMAATNRPNISIDPALRRFGRFDREVDIG IPDATGRLEILQIHTKNMKLADDVDLEQVANETHGHVGADLAALCSEAALQAI RKKMDLIDLEDETIDAEVMNSLAVTMDDFRWALSQSNPSALRETVEVPQVTW EDIGGLEDVKRELQELVQYPVEHPDKFLKFGMTPSKGVLFYGPFGCGKTLAK AIANECQANFISIKGPELLTMWFGSEANVREIFDKARQAAPCVLFFDELDSI AKARGGNIGDGGGAADRVINQILTEMDGMSTKKNVFIIGATNRPDIIIDPAILR PGRDLQLIYIPLPDEKSRVAILKANLRKSPVAKDVLDFLAKMTNGFSGADLT EICQRACKLAIRESEIEIRREERQTNPSAMEVEEDDPVEIRRDHFEEAMR FARRSVSDNIRKYEMFAQTLQQSRGFGSFRFPNGSGGAGPSQSGSGGGTGS VYTEDNDDDLGSSKGGYGLNDIFEAKIEWHE</p> |

| Protein mass | with biotin, + 226.1 | Protein pI | Extinction (M⁻¹ cm⁻¹) |
|---------------------|---------------------------------|-------------------|--|
| 91799.7 | 92025.8 | 5.12 | 42860 |

| | |
|-----------------------------------|--|
| Expression host | E. coli BL21(DE3)-R3-BirA |
| Expression medium | <p>Terrific broth (TB) with 50 µg/ml kanamycin, 50 µg/ml Streptomycin and Chloramphenicol (34 µg/ml) included in starter cultures but not in final growth.</p> <p>Biotin (10 mM stock): Add 24mg of d-biotin to 10 ml of warm (microwaved) 10mM bicine buffer (pH 8.3) and sterilize of the solution with a syringe and a 0.2 micron filter.</p> |
| Transformation and storage | Cells were transformed with the expression plasmid and selected on LB-agar plate's containing kanamycin, streptomycin and chloramphenicol. Several colonies were picked from the overnight transformation and used to inoculate cultures in LB supplemented with KAN/Cam/Str which were grown overnight and stored as glycerol (15% v/v) stock at -80°C. |

Protein information sheet

| | |
|--|--|
| Expression and <i>in vivo</i> biotinylation | <p>10ml starter cultures in LB supplemented with KAN/Cam/Str were grown overnight at 37°C with shaking at 250rpm. Next day, 10 ml was used to inoculate 1L of TB in a 2L UltraYield flask supplemented with Streptomycin 50µg/ml and Kanamycin 50µg/ml. Note: the medium need not contain chloramphenicol at this step; the pRARE2 plasmid is not lost during the growth process and also ensures a more predictable growth. Seal flask with an adhesive air porous film during cell growth. Culture was grown in a 37°C incubator with shaking at 200rpm. Cell growth was monitored until Mid-log phase was reached and cell density OD_{600nm} was approximately 2.00 – 3.00 in TB. Incubator temperature was lowered to 18°C and culture flask allowed to cool. After 30 - 45 minutes protein expression was induced by the addition of IPTG at 0.1mM final concentration. During the induction 10ml of biotin solution is added to each 1L culture on induction with IPTG. Incubation of culture was continued for 16 hours at 18°C. Next day, the cells were collected by centrifugation (15 min, 5,000 RPM in JLS-8.100 rotors). The supernatant discarded and the cell pellet transferred into a 50ml conical polypropylene tube. Cell pellet was weighed, frozen and stored at -80°C.</p> |
| Purification buffers | <p>2x lysis buffer: 100 mM HEPES, pH 7.5, 1.0 M NaCl, 10% v/v glycerol, 1mM TCEP, 10mM imidazole. Protease inhibitors (1:200 Merck cocktail IV) were added before use.</p> <p>1x lysis buffer: 50 mM HEPES, pH 7.5, 0.5 M NaCl, 10% v/v glycerol, 1 mM TCEP, 10 mM imidazole</p> <p>Wash buffer: 50 mM HEPES, pH 7.5, 0.5 M NaCl, 10% v/v glycerol, 1 mM TCEP, 40 mM imidazole</p> <p>Elution buffer: 50mM Hepes pH 7.5, 500mM NaCl, 10% Glycerol, 1mM TCEP, 100mM, 150mM, 250mM, 500 Imidazole, respectively</p> <p>GF buffer: 20 mM HEPES, pH 7.5, 0.5M NaCl, 5% v/v glycerol, 1 mM TCEP</p> <p>PEI 5%, pH 7.5: 5% polyethyleneimine (diluted from 50% stock, Sigma) and titrated to pH 7.5 with HCl.</p> |
| Purification step 1: IMAC | <p>Cell Lysis: Frozen cell pellets from -80°C were thawed by intermittently placing in a 37°C water bath, then transferred directly to ice. Cells were re-suspended in pre-chilled lysis buffer by vortex or pipetting gentle until cells were completely in suspension. 1 volume (v/w) of 2x lysis buffer was added, and the cells were thoroughly suspended. After further addition of 3 volumes (3 ml per 1g of pellet) of 1x lysis buffer the cells were broken by sonication. Lysate was then transferred to a 50ml Beckman polypropylene centrifuge tube. Cell debris and DNA from the lysate was precipitated by adding Polyethylenimine (PEI) to the tube from a 5% stock solution to a final concentration of 0.15% and incubated on ice for 15 – 30 mins for precipitation. Cell lysates were clarified by centrifugation using a JA25.50 rotor at 21,500rpm for 60 min at 4°C. The cleared supernatant was collected and filtered using a 0.45µm acrodisc syringe filter and transferred to a fresh 50ml-conical polypropylene tube.</p> <p>Ni²⁺-IDA IMAC (column 1): A 25ml Econo-Pac polypropylene column (BioRad) was packed with 5ml of Ni-Sepharose resin and equilibrated in 1x lysis buffer. The cell lysate was loaded onto the column and passed through resin at gravity flow. The unbound flow-through was collected, a</p> |

Protein information sheet

| | |
|--|---|
| <p>Step 2: Removal of His₆ tag and Ni²⁺-IDA reverse binding</p> | <p>wash step of 10 column volumes of wash buffer was included and collected as the wash fraction. Protein was eluted in four intermediate fractions with a step gradient of 5ml imidazole elution buffers; 100mM, 150mM, 250mM and 500mM imidazole, respectively, each gradient eluted fraction was collected. IMAC fractions were analysed by SDS-PAGE gel for expected molecular weight and purity, the selected eluted fractions were pooled together.</p> <p>Mass spectrometry characterisation: Identity and biotinylation efficiency was verified by mass spectrometry, using a Quadrupole-Time of flight (Q-TOF) analyser, Biotin quantitation is measured by: MW of biotinylated protein = MW of sequenced Protein + 226 Da (biotinylation modification).</p> <p>Enzymatic removal of histidine tag: Selected, pooled eluted fractions from Ni²⁺-IDA IMAC were concentrated using a 10,000 Amicon-Ultra MWCO centrifugal concentrator to approximately 10mls in volume. Cleavage of the histidine tag was performed by the addition of His₈-tagged TEV protease in a 1:20 w/w ratio to the pooled eluted protein fractions. The protein sample plus TEV was transferred to Snakeskin dialysis tubing (3,500 MWCO) and dialysed in 1L Gel Filtration buffer overnight at 4°C on a magnetic stirrer plate for digestion</p> <p>Ni²⁺-IDA Rebinding column (Column 2): To remove uncleaved protein and other contaminants, the dialysate was passed through a second IMAC column.</p> <p>1ml of packed Ni-Sepharose resin was added to a 25ml Econo-Pac polypropylene column (BioRad), washed in 3 CV of distilled deionized water then equilibrated in 5 column volumes of gel filtration buffer. The cleaved, dialyzed protein sample was loaded and passed through the column by gravity flow. The flow through was collected and combined with two 2 CV washes in GF buffer gel filtration buffer. For completeness, the column was finally eluted with a 2 CV wash of 500mM imidazole elution buffer. All fractions were analysed for verification of cleaved protein and purity by SDS-PAGE gel and mass spectrometry characterisation.</p> |
| <p>Step 3: SEC: Gel filtration</p> | <p>Size Exclusion Chromatography (SEC) (Column 3): Verified protein from rebinding column 2 was concentrated to required volume for gel filtration (<5 ml), using an Amicon-Ultra 30,000 MWCO centrifugal concentrator. The GF column (Superdex S200 16/60 HiLoad, 120ml on an ÄKTA Express chromatography system) was equilibrated in gel filtration buffer. The concentrated protein sample was loaded onto the GF column for further purification; protein was fractionated on the column in GF buffer at 1ml/min and 1ml fractions were collected at the A280 peaks.</p> <p>Peak fractions were analysed by SDS-PAGE gel, some impurities and contaminating bands were observed on the gel analysis. Peak GF fractions containing least impurities were collected, pooled and concentrated for an additional second GF purification step, attempting to reduce contaminants and increase purity of protein to >90%. Protein purity was improved after additional GF step with protein purity increasing to 90%. Selected peak GF fractions containing high protein concentration and purity were collected and pooled for the protein final sample.</p> <p>Protein concentration mg/ml and storage: Protein was concentrated using a 30,000 Amicon-Ultra MWCO centrifugal concentrator and final</p> |


Protein information sheet



| | |
|--------------------------------|--|
| Step 4: Protein storage | protein concentration was determined from the absorbance at 280nm, measured using a nanodrop spectrometer (Nanodrop ND-1000, NanoDrop Technologies). Final protein concentration was measured at 1.3mg/ml with a final volume of 1.6mls. Protein was transferred into 200 μ l aliquots, flash-frozen in liquid nitrogen and stored at -80°C. The final storage buffer was 20mM Hepes, pH 7.5, 500mM NaCl, 5% v/v glycerol, 1 mM TCEP |
|--------------------------------|--|


Gel Filtration Chromatogram (Superdex S200)

VCPA

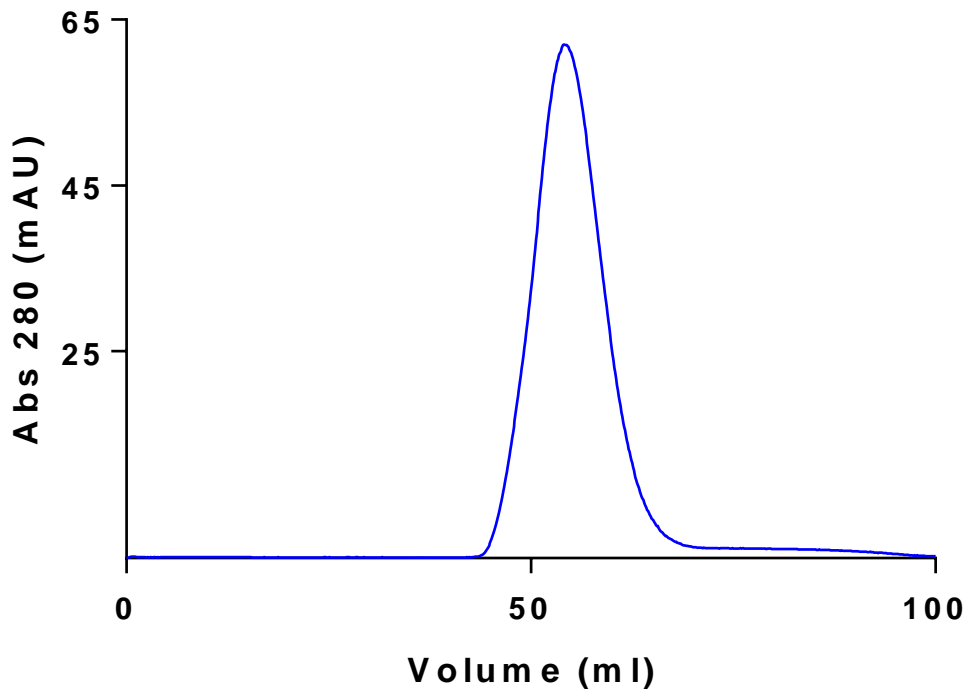

VCPA-c000, Maximum peak elution time was 54.241ml with absorbance at 280nm of 61.836 mAU. SEC analysis validates protein is a dimer.

Final Sample SDS-PAGE

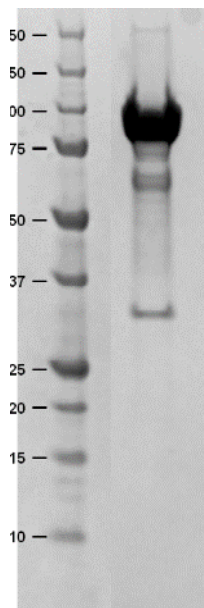

| Protein | Expected mass + Biotinylation (Da) | Final concentration (mg/ml) |
|-----------|------------------------------------|-----------------------------|
| VCPA-c000 | 92025.8 | 1.3 |

Protein information sheet



Mass Spectrometry characterisation

| Protein | Vector | Expected mass cleaved + Biotinylation (Da) | Observed mass | Comments |
|-----------|-----------|---|------------------|-------------------|
| VCPA-c000 | pNIC-Bio3 | 92025.8 | 92026 | 100% Biotinylated |

